

TRAPSYSTEM® SET

SET PER BIOPSIA OSTEOMIDOLLARE CON DISPOSITIVO TRAPSYSTEM®
BONE MARROW BIOPSY SET WITH TRAPSYSTEM® DEVICE

TRAPSYSTEM® SET è l'ago per biopsia osteomidollare che ha rivoluzionato la biopsia ossea eliminando la manovra di lussazione.

Caratteristiche del dispositivo:


- ago tipo Jamshidi rastremato, centimetrato, con punta a bocca di pesce con impugnatura ergonomica con connettore luer lock che consente di distribuire la forza su tutto il palmo della mano e non localmente
- mandrino estraibile, perfettamente accoppiato all'ago, con punta piramidale a tre facce
- dispositivo Trapsystem, privo di punta, disegnato per scorrere all'interno dell'ago Jamshidi, caratterizzato da una scanalatura che consente di imprigionare il frustolo, eliminando la manovra di lussazione
- filo di estrazione a punta smussa, centimetrato, indicato per la valutazione dell'entità del prelievo prima della biopsia e per il recupero del campione
- dispositivo guida di sicurezza indicato per ricoprire la punta dell'ago ed evitare pericoli per l'operatore
- tappo luer lock indicato per chiudere l'ago evitando la fuoriuscita di sangue

TRAPSYSTEM® SET is the first bone marrow biopsy needle that avoids the deflection technique.

Device Features:

- Jamshidi needle with tapered "fish mouth" tip, depth marks and ergonomic handle with luer lock connector (it ensures safe grip and manoeuvrability and force distribution on the palm of the hand)
- removable stylet with pyramidal tip
- Trapsystem device with tip designed to trap the sample into the Jamshidi needle, avoiding the deflection technique
- extraction wire with blunt tip intended for sample size evaluation and sample retrieval.
- safety guide intended for needle tip covering to avoid damages for personnel
- luer lock cap

CODE	DIAMETER	LENGTH	BOX
TRAPJ 0810	8G	100mm	10
TRAPJ 0815	8G	150mm	10
TRAPJ 1110	11G	100mm	10
TRAPJ 1115	11G	150mm	10
TRAPJ 1306	13G	60mm	10
TRAPJ 1310	13G	100mm	10

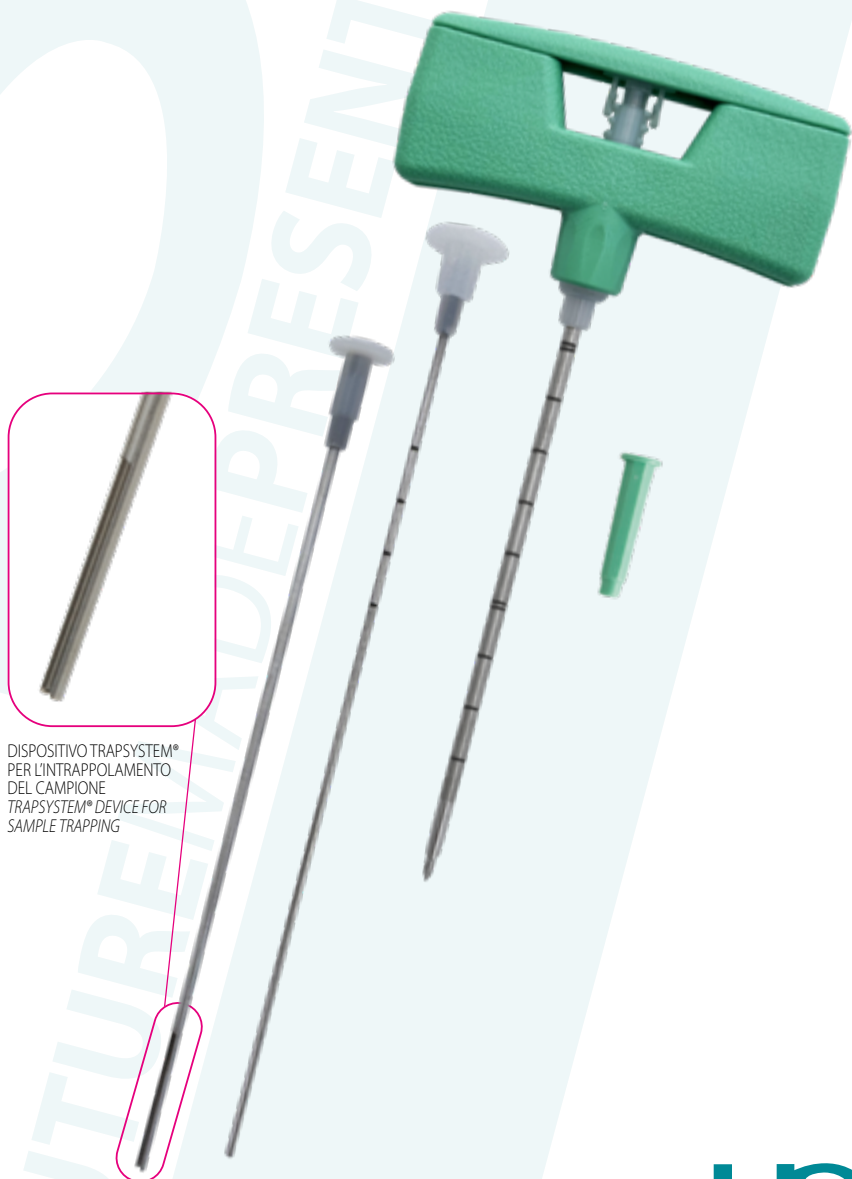

DISPOSITIVO TRAPSYSTEM®
PER L'INTRAPPOLAMENTO
DEL CAMPIONE
TRAPSYSTEM® DEVICE FOR
SAMPLE TRAPPING